# Julabo

# DYNEO DD-601F Refrigerated / heating circulator

Refrigerated circulators of the DYNEO series distinguish themselves with a great price-performance ratio. The instruments offer high heating/cooling capacities for short heat-up and cool-down times. The refrigerated circulators work precisely and reliably even at higher ambient temperatures up to +40 °C. Either in basic research, in material testing or in technical systems – the DYNEO refrigerated circulators offer functional solutions for every requirement and budget.


### Optional analog and digital interface

DYNEO circulators can optionally be equipped with analogue or digital interfaces. To request the options, order number must be extended with .d for the digital and .a for the analog interface (9XXX XXXX.A / 9XXX XXX.D)




### **Product features**

- Optimized cooling coil design saves space in the bath tank
- powerful and infinitely adjustable pressure pump
- Flow rate 27 l/min, pressure 0.7 bar
- · easy switching between internal and external circulation
- large color TFT display, multilingual interface
- central rotary knob (controller) simplifies operation
- Integrated programmer
- Integrated external Pt100 connection
- USB connection
- RS232 interface or analog interfaces (optional)
- integrated drain makes emptying liquid easy and safe.
- Bath cover included with delivery
- removable ventilation grid
- Powerful cooling machines
- For internal and external applications
- Integrated pump connection M16×1
- For internal and external applications

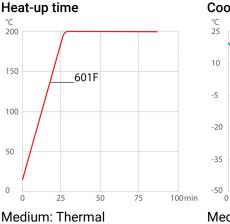

### Cool-down time

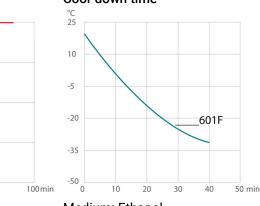

Medium: Ethanol





### Performance values

230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)		
Heating capacity kW	2	
Viscosity max. cSt	50	
Pump capacity flow rate l/min	8 27	
Pump capacity flow pressure bar	0.1 0.7	
Power consumption A	14	

Order No.			9021705.33		
Cooling capacity (Ethanol	)				
°C	20	0	-10	-20	-30
kW <sup>1</sup>	0.6	0.44	0.27	0.16	0.04
Refrigerant stage 1					
Refrigerant	R449A				
Filling weight g	150				
Global Warming Potential R449A	for 1397				
Carbon dioxide equivalent	tt 0.2095	5			

<sup>1</sup> Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

### **Technical data**

Available voltage ve	ersions		Cooling	
Order No.	9 021 705		Cooling of compressor	1-stage Air
Available voltage versions:				
9021705.01	100V/50-60Hz (N (R452A)	lema N5-15 Plug)		
9021705.02	5.02 115V/60Hz (Nema N5-15 Plug) (R449)			
9021705.33	200-230V/50-60H 7/4 Plug Type F)	Hz (Schuko Plug - CEE (R449A)		
9021705.04	200-230V/50-60 BS1363A) (R449			
9021705.05	200-230V/50-60H 1011) (R449A)	Hz (CH Plug Type SEV		
9021705.33.chn	200-230V/50-60H	Hz (CN Plug) (R449A)		
Bath			Other	
Bath tank		Stainless steel	Classification	Classification III (FL)
Bath cover		integrated	Pump function	Pressure Pump
Usable bath opening cm (W x L / D) 22 x 15 / 20		22 x 15 / 20	Pump type	Immersion Pump
			User Interface Language	Chinese, English, French, German, Italian, Japanese, Korean, Portuguese, Russian, Spanish
Electronics			Dimensions and volumes	
			Weight kg	38.2



8/12 mm 33 x 47 x 74 8 ... 10 M16x1 male

	Interfaces	Alarm output optional,	Barbed fittings inner diameter mm
		REG/EPROG optional, RS232 optional, Standby-Input	Dimensions cm ( $W \times L \times H$ )
			Filling volume I
	optio	optional, USB	Pump connections
	External pt100 sensor connection	integrated	
	Integrated programmer	8x60 steps	
	Temperature control	PID3	
	Absolute temperature calibration	3 Point Calibration	
	Temperature display	3.5" TFT Display	
	Temperature setting	Shaft Encoder	
	Electronic Timer h:min	00:00 99:59	

Temperature values	
Setting the resolution of the temperature display °C	0.01
Working temperature range °C	-35 +200
Temperature stability °C	±0.01
Ambient temperature °C	+5 +40
Setting the resolution of the temperature display °C	0.01

### Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male)

### All Benefits



### Handle with ease.

Makes day-to-day work easy. Comfortably move your JULABO Circulator around by using the ergonomic handles (front and rear).



### Wide range.

Refrigerated and heating circulator in various combinations, circulator in various sizes. Maximum flexibility through a large selection of accessories.



Information. Everything clear. Information in plain text on a large color screen.



**Turn. Push. Go.** Easy operation of all parameters using the central controller.



Powerful. Adjustable. Strong pressure pump, continuously adjustable.



# Highly precise

PID Temperature control with drift compensation and adjustable control parameters, temperature stability  $\pm 0.01...\pm 0.02$  °C



**Brilliance. In color.** Large color display with vivid luminance is easy to read, even from a large distance.



Multi-lingual. Operation in multiple languages.



**Programmer. Integrated.** The integrated internal programmer makes it possible to automatically run temperature time profiles.



USB. Remote control made easy using the integrated USB interface.





**RS232.** Connection using the optional RS232 interface.



Temperature. Under control. External Pt100 sensor connection for precise measurement and control directly in the external application.



**Process stability.** Early warning - visual and acoustic - of critical states increases process stability.



ATC3. Calibration.

'Absolute Temperature Calibration' for compensating a physically caused temperature difference, 3-point calibration.



# Connection. Easy.

Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



# More bath.

Designed for more comfort. Thanks to the recessed cooling coil, the internal bath provides more space.



## Solid. Minimized energy loss through high-quality insulation.



### Condensation protection.

Superb design solution. Integrated ventilation directs air over the bath lid and minimizes condensation.



**Green technology.** Development consistently applied environmentally

friendly materials and technologies.



Quick start. Individual JULABO consultation and comprehensive manuals at your disposal.



Analog I/O. Analog interfaces for integration into process control systems (optional).



Fill level. Monitored. Fill level indicator on the display for heat-transfer liquid.



# Process. Under control.

Full regulation of the dynamics control, access to all important control parameters for individual process optimization.



### Stable. Mobile.

Rubber feet keep JULABO Circulators standing firm. Larger and more powerful units also have integrated rollers for easy handling.



### 100 % Cooling capacity

'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher temperatures



### Space saving. Free up space.

Place your JULABO Circulator right next to an application, another unit, or wall. That saves space. This is made possible by eliminating vents and connections on the sides.



### Tidy.

The special drain tap for easy draining of bath fluids without tools.



### 100% Checked.

100% testing. 100% quality. Each JULABO Circulator undergoes thorough quality testing before leaving the factory.



JULABO. Quality. Highest standards of quality for a long product life.



Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.





#### Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.



#### Quick support

If an error occurs, the integrated Black-Box function permits fast diagnosis by the JULABO service team