

## CORIO CD-1201F Refrigerated / heating circulator

Refrigerated / heating circulators of the CORIO CD series distinguish themselves with a great price-to-performance ratio. Durable and dependable, they are suitable for a wide range of standard and routine work in both laboratories and industrial environments.

The CORIO CD-1201F refrigerated circulator is exceptionally versatile with a bath opening of 35 x 41 cm (width x depth) and a bath depth of 30 cm. The resulting bath volume of 50 liters makes it possible to control the temperature of large objects, or many smaller objects, directly in the bath.


This cooling machine uses natural, environmentally-friendly refrigerant and was developed with a focus on energy efficiency. This means significant savings of up to 70% with respect to the operating costs for numerous applications, which also means faster amortization of the procurement cost. At the same time, the lower energy consumption positively contributes to climate protection.



### Product features

- Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- USB connection
- Integrated pump connection M16×1
- easy switching between internal and external circulation
- Bath lid and drain tap included
- removable ventilation grid
- Refrigeration unit without side vents
- Class III (FL) according to DIN 12876-1
- Space-saving cooling coil design yields more usable space in the bath tank

### Pump capacity

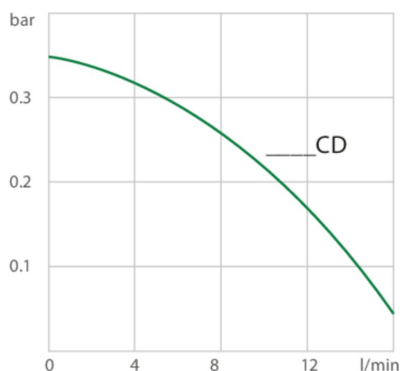

Medium: Water

## Performance values

### 230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)

Heating capacity kW	2
Viscosity max. cSt	50
Pump capacity flow pressure l/min	15
Pump capacity pressure bar	0.35
Power consumption A	11

Order No. 9012718.N1.03

#### Cooling capacity (Ethanol)

°C	20	10	0	-10	-20	-30	-40
kW <sup>1</sup>	1.25	1.2	1.1	0.9	0.63	0.38	0.21

#### Note about natural refrigerants:

Temperature control units using natural refrigerants are often subject to regulatory requirements regarding the installation site, operation, transport or disposal of the units. If you have any questions, we will be happy to advise you.

#### Refrigerant stage 1

Refrigerant	R1270
Filling weight g	90
Global Warming Potential for R1270	2
Carbon dioxide equivalent t	0.00018

<sup>1</sup> Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

## Technical data

### Available voltage versions

Order No. 9 012 718

#### Available voltage versions:

9012718.N1.03	230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F) (R1270)
9012718.N1.03.chn	230V/50Hz (CN Plug) (R1270)
9012718.N1.04	230V/50Hz (UK Plug Type BS1363A) (R1270)
9012718.N1.05	230V/50Hz (CH Plug Type SEV 1011) (R1270)
9012718.N1.13	230V/60Hz (Nema N6-20 Plug) (R1270)
9012718.N1.01	100V/50-60Hz (Nema N6-20 Plug) (R1270)
9012718.N1.02	115V/60Hz (Nema N6-20 Plug) (R1270)

### Cooling

Cooling of compressor 1-stage Air

### Bath

Bath tank	Stainless steel
Bath cover	integrated
Usable bath opening cm (W x L / D)	35 x 40 / 30

### Other

Classification	Classification III (FL)
IP Code	IP 21
Pump function	Pressure Pump
Pump type	Immersion Pump

Electronics	
Interfaces	USB
Temperature control	PID1
Absolute temperature calibration	1 Point Calibration
Temperature display	LED
Temperature setting	Keypad

Dimensions and volumes	
Weight kg	68
Barbed fittings inner diameter mm	8/12 mm
Dimensions cm (W x L x H)	45 x 64 x 95
Filling volume l	48 ... 56
Pump connections	M16x1 male

Temperature values	
Setting the resolution of the temperature display °C	0.1
Working temperature range °C	-40 ... +100
Temperature stability °C	±0.03
Ambient temperature °C	+5 ... +40
Temperature display resolution °C	0.01 ... 0.1

**Included in delivery**  
 2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male)

### All Benefits



**Brilliant.**  
 Very bright display makes it easy to read even from a distance.



**Switch on. And off you go.**  
 Intelligent operating concept. Ready for operation with just a few quick and easy steps.



**Timer. Integrated.**  
 CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes.



**Everything at the front.**  
 All operating controls and safety functions are accessed easily and comfortably from the front.



**Exact.**  
 You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO precise and perfect.



**Locked in.**  
 The lockable power plug guarantees safe connection. More process safety.



**Connection. Easy.**  
 Inclined pump connections (M16x1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



**Early warning system for low liquid level**  
 Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



**Connectivity.**  
 Remote control made easy. CORIO CD circulators feature a USB connection.